

CyanoKit Training Materials

Training Kit includes (one of each): Glass medication vial, Transfer Spike, IV Solution Set, 100 mL NS (will use a 250 mL bag for the actual Cyanokit), alcohol swab

1. **Open the top of the glass vial to expose the rubber stopper and clean with an alcohol swab**
2. **Open the transfer spike**
3. **Remove the rubber stopper on the end of the NS bag (the rubber stopper is blue in the training kit)**




4. **Uncap the slim end of the transfer spike and insert this into the exposed port of the NS bag**



5. **Uncap the other end of the transfer spike and insert into the middle of the rubber stopper in the glass vial**
 - **In the training materials, the glass vial is not a vacuum vial, therefore the solution will need to be squeezed into the glass vial**
 - **In the CyanoKit, the glass vial is a vacuum vial and the solution will readily be pulled into the glass vial**

6. In the CyanoKit, the glass vial should be filled to the line

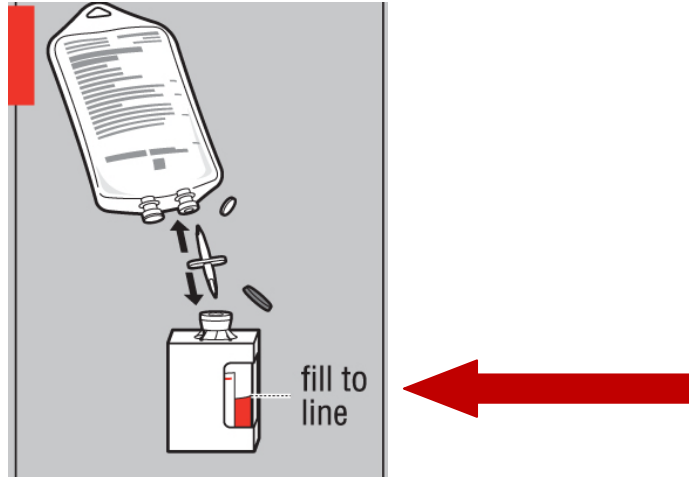

7. Remove the transfer spike and swirl the contents until reconstituted
8. Uncap and insert the IV Solution spike into the same hole (that was made by the transfer spike) in the rubber stopper of the glass vial
9. Invert and prime the tubing. If your IV Solution set has a vent, open the vent.

